Check for updates

OPEN ACCESS

EDITED BY Anna B. Nilsson, Karolinska Institutet (KI), Sweden

REVIEWED BY

Gianluca Vergine, Ospedale Infermi di Rimini, Italy Samuel Arthur Rhedin, Karolinska Institutet (KI), Sweden

*CORRESPONDENCE

Alexis Rybak alexis.rybak@activ-france.fr Corinne Levy corinne.levy@activ-france.fr

SPECIALTY SECTION

This article was submitted to Pediatric Infectious Diseases, a section of the journal Frontiers in Pediatrics

RECEIVED 28 June 2022 ACCEPTED 29 August 2022 PUBLISHED 21 September 2022

CITATION

Cohen R, Aupiais C, Filleron A, Cahn-Sellem F, Romain O, Béchet S, Auvrignon A, Batard C, Virey B, Jung C, Rybak A and Levy C (2022) Diagnostic accuracy of SARS-CoV-2 rapid antigen test from self-collected anterior nasal swabs in children compared to rapid antigen test and RT-PCR from nasopharyngeal swabs collected by healthcare workers: A multicentric prospective study. *Front. Pediatr.* 10:980549. doi: 10.3389/fped.2022.980549

COPYRIGHT

© 2022 Cohen, Aupiais, Filleron, Cahn-Sellem, Romain, Béchet, Auvrignon, Batard, Virey, Jung, Rybak and Levy. This is an open-access article distributed under the terms of the Creative Commons Attribution License (CC BY). The use, distribution or reproduction in other forums is permitted, provided the original author(s) and the copyright owner(s) are credited and that the original publication in this journal is cited, in accordance with accepted academic practice. No use, distribution or reproduction is permitted which does not comply with these terms.

Diagnostic accuracy of SARS-CoV-2 rapid antigen test from self-collected anterior nasal swabs in children compared to rapid antigen test and RT-PCR from nasopharyngeal swabs collected by healthcare workers: A multicentric prospective study

Robert Cohen^{1,2,3,4,5}, Camille Aupiais^{6,7}, Anne Filleron⁸, Fabienne Cahn-Sellem⁴, Olivier Romain^{4,9}, Stéphane Béchet^{1,4}, Anne Auvrignon⁴, Christophe Batard⁴, Brigitte Virey⁴, Camille Jung^{2,3}, Alexis Rybak^{1,4,6,10*} and Corinne Levy^{1,2,3,4,5*}

¹Association Clinique et Thérapeutique Infantile du Val-de-Marne (ACTIV), Créteil, France, ²Clinical Research Center, Centre Hospitalier Intercommunal de Créteil, Créteil, France, ³Université Paris Est, IMRB-GRC GEMINI, Créteil, France, ⁴Association Française de Pédiatrie Ambulatoire (AFPA), Orléans, France, ⁵Groupe de Pathologie Infectieuse Pédiatrique (GPIP), Créteil, France, ⁶Assistance Publique–Hôpitaux de Paris, Clinical Epidemiology Unit, Robert Debré University Hospital, ECEVE INSERM UMR 1123, Université de Paris, Paris, France, ⁷Assistance Publique–Hôpitaux de Paris, Pediatric Department, Jean Verdier University Hospital, Université Sorbonne Paris Nord, Bondy, France, ⁸Department of Pediatrics, CHU Nîmes, University of Montpellier, Nîmes, France, ⁹Assistance Publique–Hôpitaux de Paris, Neonatalogy Department, Antoine Béclère Universital, Université Paris, Pediatric Emergency Department, Robert Debré University Hospital, Université de Paris, Pediatric

Testing for SARS-CoV-2 is central to COVID-19 management. Rapid antigen test from self-collected anterior nasal swabs (SCANS-RAT) are often used in children but their performance have not been assessed in real-life. We aimed to compare this testing method to the two methods usually used: reverse transcription polymerase chain reaction from nasopharyngeal swabs collected by healthcare workers (HCW-PCR) and rapid antigen test from nasopharyngeal swabs collected by healthcare workers (HCW-PCR), and rapid antigen test from nasopharyngeal swabs collected by healthcare workers (HCW-RAT), estimating the accuracy and acceptance, in a pediatric real-life study. From September 2021 to January 2022, we performed a manufacturer-independent cross-sectional, prospective, multicenter study involving 74 pediatric ambulatory centers and 5 emergency units throughout France. Children \geq 6 months to 15 years old with suggestive symptoms of COVID-19 or children in contact with a COVID-19–positive patient were prospectively enrolled. We included 836 children (median 4 years), 774 (92.6%)

were symptomatic. The comparators were HCW-PCR for 267 children, and HCW-RAT for 593 children. The sensitivity of the SCANS-RAT test compared to HCW-RAT was 91.3% (95%CI 82.8; 96.4). Sensitivity was 70.4% (95%CI 59.2; 80.0) compared to all HCW-PCR and 84.6% (95%CI 71.9; 93.1) when considering cycle threshold <33. The specificity was always >97%. Among children aged \geq 6 years, 90.9% of SCANS-RAT were self-collected without adult intervention. On appreciation rating (from 1, very pleasant, to 10, very unpleasant), 77.9% of children chose a score \leq 3. SCANS-RAT have good sensitivity and specificity and are well accepted by children. A repeated screening strategy using these tests can play a major role in controlling the pandemic.

KEYWORDS

COVID-19, ambulatory setting, test, pediatric, self-test

Introduction

Following the successive COVID-19 waves due to several SARS-CoV-2 variants, in many countries, healthcare authorities implemented non-pharmaceutical interventions and large-scale testing strategies (1, 2). Two methods were mainly used without distinction in France: reverse transcription polymerase chain reaction from nasopharyngeal swabs collected by healthcare workers (HCW-PCR), and rapid antigen test from nasopharyngeal swabs collected by healthcare workers (HCW-PCR), in addition to immunization programs (3). In 2021, 48.8% of the 168 million of tests recorded in the French national database were HCW-RAT (4). While HCW-RAT has lower analytical sensitivity than HCW-PCR, this method is highly specific, inexpensive, and provides results in minutes.

Testing for SARS-CoV-2 is central to COVID-19 management and essential to detect people who are likely infectious, helping to implement control measures (5). For SARS-CoV-2 testing and screening, especially in children, rapid antigen test from self-collected anterior nasal swabs (SCANS-RAT) could be a useful tool (5). The duration of the pandemic and the frequency with which testing must be done to limit infectiousness, particularly in schools, means that the less invasive, less painful and less unpleasant tests should be used, to avoid poor acceptance by children and families.

The diagnostic accuracy of HCW-RAT for diagnosing SARS-CoV-2 infection in children has been assessed in several studies and was the subject of a meta-analysis (6). No test included fully satisfied the performance requirements recommended by the World Health Organization, and the diagnostic accuracy of the HCW-RAT under real-life conditions varied broadly (6). A recent French study in an emergency department found good sensitivity of the HCW-RAT in real life for symptomatic children, and when focused on high viral load, the sensitivity was excellent (7). Nasopharyngeal swabbing compared to other upper-respiratory sampling methods, including oropharyngeal swab, appeared to be superior in a pediatric study finding a significantly higher positivity rate and a significantly higher viral load on nasal samples (8). In adults, the diagnostic accuracy of SCANS-RAT was assessed in several studies, but relatively few patients were enrolled (9-12). Millions of SCANS-RAT are used each day worldwide, but to our knowledge, no study has assessed their performance in real-life in children.

This study compared SCANS-RAT to HCW-PCR and HCW-RAT, estimating the accuracy and acceptance, in a pediatric real-life study.

Methods

From September 10, 2021, to January 29, 2022, the Association Clinique et Thérapeutique Infantile du Val de Marne (ACTIV) network conducted a manufacturer-independent cross-sectional, prospective, multicenter study involving 74 pediatric ambulatory centers (see the Acknowledgments section) and 5 emergency units (Jean-Verdier hospital in Bondy, intercommunal hospital of Créteil, Princess Grace hospital in Monaco, Carémeau hospital in Nîmes, and Versailles hospital in Le Chesnay) throughout France. Children \geq 6 months to 15 years old with suggestive symptoms of COVID-19 or children in contact with a COVID-19–positive patient were prospectively enrolled.

Abbreviations: Ct, cycle threshold; HCW-RAT, rapid antigen test from nasopharyngeal swabs collected by healthcare workers; HCW-PCR, reverse transcription polymerase chain reaction from nasopharyngeal swabs collected by healthcare workers; SCANS-RAT, rapid antigen test from self-collected anterior nasal swabs.

Healthcare workers collected nasopharyngeal swabs to perform a rapid antigen test during the medical visit. Reverse transcription polymerase chain reaction were collected either during the medical visit or with a medical prescription in a laboratory. Swabs were performed as recommended in international guidelines (13). Ambulatory and hospital virology laboratories analyzed the HCW-PCR specimens according to the French National Reference Center recommendations (14). At the same time, SCANS-RAT was offered to children in the pediatrician office or emergency department. After oral instructions from an adult (parents or pediatricians), children self-collected the nasal specimen from both nares. Adults could help to perform the test when the children were not able to perform the swabbing alone. The test used was COVID-VIRO ALL IN^(K) (AAZ-LMB, Boulogne Billancourt, France) which has a short soft sponge sampling part (1.5 cm). Recommended sampling duration was 30 s (15 s per nostril) (15). Children were asked to rate the SCANS-RAT from 1, very pleasant, to 10, very unpleasant.

After informing the parents of the participating children about the study, an electronic case report form in a secure database was prospectively completed by the pediatrician. Any child or parent had the right to object to the data collection for this study.

The diagnosis accuracy of the SCANS-RAT was compared with that of HCW-PCR and/or HCW-RAT. According to the spread of different variants in France, we defined 2 periods: period 1, when the Delta variant was predominant and Omicron not yet or poorly isolated in France (from September 10, 2021 to December 19, 2021), and period 2, when the Omicron variant was spreading and became predominant (i.e., > 50%, from December, 20, 2021 to January 29, 2022) (16). We performed an ad-hoc subgroup analysis on children who had a HCW-PCR with Ct <33 and with Ct <30. Data were entered by using an electronic case report form (PHP/MySQL) and were analyzed by using Stata/SE v15 (StataCorp, College Station, TX, USA). Quantitative data were compared by Student *t* test and categorical data by chi-squared or Fisher exact test.

Results

Among the 836 patients with a SCANS-RAT (median 4.0 years, interquartile range 2-7), 774 (92.6%) were symptomatic. In addition to a SCANS-RAT, 263 children had a HCW-PCR, 589 children a HCW-RAT, and 4 children both tests. Patients characteristics are detailed in Table 1. The prevalence of SARS-CoV-2 infection was 18.4% (154/836) (95%CI 15.8; 21.2) during the whole study period: 10.1% (62/617) (95%CI 7.8; 12.7) in period 1 (Delta wave) and 42.0% (92/219) (95%CI 35.4–48.8) in period 2 (Omicron wave).

The cycle threshold (Ct) was available for 75.3% (61/81) of positive HCW-PCR. The Table 2 shows the performance of the SCANS-RAT compared to HCW-PCR and HCW-RAT. The overall sensitivity of the SCANS-RAT compared to any positive HCW-PCR results was 70.4% (95%CI 59.2; 80.0). False negative SCANS-RAT results compared to HCW-PCR (9.0%, n = 24/267) corresponded mainly to HCW-PCR tests with Ct \geq 30 (Figure 1). Thus, sensitivity was 84.6% (95%CI 71.9; 93.1) and 93.6% (95%CI 82.5; 98.7) when considering only HCW-PCR with Ct <33 and with Ct <30, respectively. The specificity was always high (from 97.4 to 97.8%). The median delay between the SCANS-RAT and the HCW-PCR was 0 day (interquartile range 0–0). This delay was similar between positive and negative HCW-PCR (p = 0.50).

The sensitivity of the SCANS-RAT compared to HCW-RAT results was 91.3% (95%CI 82.8; 96.4) while specificity was 99.6% (95%CI 98.6; 100). Of note, sensitivity was comparable between period 1 (85.2%, 95%CI 66.3; 95.8) and period 2 (94.3%, 95%CI 84.3; 98.8). False negative SCANS-RAT results compared to the HCW-RAT accounted for 7/593 cases.

For children from ≥ 2 to <6 years old (n = 308), the SCANS-RAT was obtained by the children themselves alone in 47.4% (146/308) of cases, while for 52.6% (162/308) of cases, an adult, mostly the accompanying adult, had to help. From age ≥ 6 years, the self-collecting swab was easily performed by the children themselves in 90.9% (179/197) of cases. On appreciation rating (from 1, very pleasant, to 10, very unpleasant), 77.9% of children chose a score ≤ 3 .

Discussion

To our knowledge, this is the first large prospective study assessing in real-life the diagnostic accuracy of an SCANS-RAT in a pediatric population. Studies assessing SCANS-RAT diagnosis accuracy have involved only adults (9-12). As compared with the HCW-RAT, for which a positive test result is often associated with live viral culture (17), the diagnostic accuracy of the SCANS-RAT is similar in both sensitivity and specificity. If we consider all positive tests independent of Ct number, as compared with HCW-PCR, the SCANS-RAT had excellent specificity but relatively moderate sensitivity under the minimum performance requirements as recommended by the World Health Organization (6). However, HCW-PCR have been reported to remain positive up to 5 weeks after infection while live virus is usually isolable only during the first week (5). Thus, HCW-PCR without Ct may not be the gold-standard to detect contagious patients. There is a continuous relation between Ct and viral culture with a 33% reduction of the odds of live viral culture for 1-unit increase in Ct (18). Several thresholds have been proposed (18). In France, tests with $Ct \leq 33$ are reported as "positive" whereas tests with Ct > 33 are reported as

TABLE 1 Characteristics of patients included in the study (N = 836).

Test performed	HCW-RAT $N = 593$	HCW-PCR $N = 267$	SCANS-RAT $N = 836$
Age in years, median (IQR)	4 (2-7)	4 (2-8)	4 (2-7)
Sexe (male), <i>n</i> (%)	321 (54.1%)	148 (55.4%)	460 (55.0%)
Positive test, <i>n</i> (%)	80 (13.5%)	81 (30.3%)	131 (15.7%)
Contact with confirmed COVID-19 positive case and asymptomatic, $n(\%)$	3 (9.1%)	3 (12.5%)	6 (10.9%)
Contact with confirmed COVID-19 positive case and symptomatic, n (%)	42 (40.4%)	49 (52.1%)	76 (40.0%)
No contact with confirmed COVID-19 positive case and symptomatic, $n(\%)$	34 (7.5%)	28 (19.3%)	48 (8.2%)

4 children had both test, data on contact with confirmed COVID-19 case was missing for 7 children.

TABLE 2 Performance of the rapid antigen test from self-collected anterior nasal swabs (SCANS-RAT) test compared to reverse transcription polymerase chain reaction from nasopharyngeal swabs (HCW-PCR), and rapid antigen test from nasopharyngeal swabs collected by healthcare workers (HCW-RAT) (N = 836).

SCANS-RAT test compared to	Sensitivity	Specificity	PPV	NPV	LR+	LR-	Accuracy
HCW-RAT	91.3	99.6	97.3	98.6	234.0	0.09	98.5
(80/593)	(82.8; 96.4)	(98.6; 100)	(90.7; 99.7)	(97.2; 99.5)	(58.6; 935)	(0.04; 0.18)	(97.1; 99.3)
All HCW-PCR	70.4	97.8	93.4	88.3	32.7	0.30	89.5
(81/267)	(59.2; 80.0)	(94.6; 99.4)	(84.1; 98.2)	(83.2; 92.4)	(12.3; 87.2)	(0.22; 0.42)	(85.2; 92.9)
HCW-PCR with Ct	84.6	97.4	89.8	96.0	33.0	0.16	94.7
<33 (52/267)	(71.9; 93.1)	(94.1; 99.2)	(77.8; 96.6)	(92.2; 98.2)	(13.8; 79.0)	(0.08; 0.30)	(91.2; 97.2)
HCW-PCR with Ct	93.6	97.5	89.8	98.5	37.4	0.07	96.8
<30 (47/267)	(82.5; 98.7)	(94.3; 99.2)	(77.8; 96.6)	(95.6; 99.7)	(15.7; 89.3)	(0.02; 0.20)	(93.7; 98.6)


Data are values (95% confidence interval).

Ct, cycle threshold; PPV, positive predictive value; NPV, negative predictive value; LR+, positive likelihood ratio; LR-, negative likelihood ratio.

"weak positive" (19). In our study, if we consider only patients with Ct < 33, which suggests a high viral load, the sensitivity was good [84.6% (95%CI 71.9; 93.1)]. Viral load is an important determinant of disease transmission, which is a critical parameter for implementing control measures and disease modeling (20, 21). The purpose of the SCANS-RAT is more to detect the most infectious patients than to accurately diagnose COVID-19.

Rapid antigen tests have multiple advantages: suitability, speed of the results and cost. Furthermore, tests from anterior nasal swabs are suitable for repeated tests in children. Indeed, during the successive epidemic waves, children had to undergo many tests, sometimes for a short period of few weeks, and good acceptability is a crucial goal: lower sensitivity of individual tests can be compensated for by frequency of testing and wider dissemination of tests. Because children show substantially reduced mortality from COVID-19, entry screening into schools might require greater compromise that balances resources and sensitivity to testing as many individuals as possible. Because of a high specificity, the risk of a false positive test due to repeated SCANS-RAT is low. The use of tests from self-collected anterior nasal swabs and not from nasopharyngeal swabs collected by healthcare workers is the first step to succeed in a largescale testing strategy allowing for widespread school opening. Repeated use of the SCANS-RAT can contribute to a wider opening of schools with expected benefits for the mental and physical health of children (22). In the United States, many schools offered free COVID-19 tests (23). In France, in early January 2022, with the Omicron wave, the testing strategy for children at school was difficult to perform: 3 tests in 5 days (24). Indeed, in this context, even if it means losing slightly sensitivity, it appears crucial to have a very good acceptance of the tests in children, allowing a wide use within families without healthcare workers support. Of note, the sensitivity of the SCANS-RAT compared to HCW-PCR and HCW-RAT did not change significantly during the delta and Omicron periods. Similar results were recently reported in a study mainly in adult population (25).

Our study has some limitations. First, we did not use centralized reverse transcription polymerase chain reaction performed by centralized high-complexity laboratories and the Ct number was available for only three-quarters of SARS-CoV-2–positive patients. However, this limitation is also a strength of our real-life study: we compared the SCANS-RAT with the

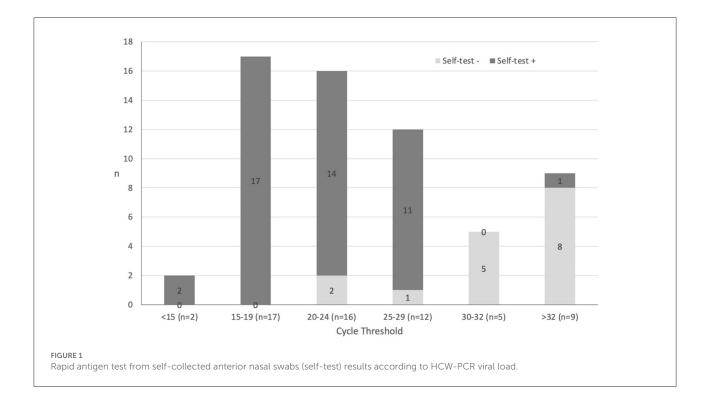

methods used in real life. Second, most children in our study were symptomatic (92.7%), and we cannot extrapolate our results for screening in asymptomatic children. However, the accuracy of the HCW-RAT and contagiousness are believed to be mainly driven by the viral load and SCANS-RAT used in our study has as good sensitivity as HCW-PCR with low Ct (5). Third, for some children, HCW-RAT and HCW-PCR were not performed the same day. However, the majority of HCW-PCR were performed in symptomatic children and in the first few days after the symptom onset. This corresponds to a period where children have high viral loads with Ct <30 (26).

In conclusion, the anterior nasal self-collected test used in this study seems reliable and suitable, allowing to detect infectious children. A repeated screening strategy using SCANS-RAT can play a major role in controlling the pandemic.

Data availability statement

The data that support the findings of this study are available from the corresponding author upon reasonable request.

Ethics statement

The study protocol was approved by an Ethics Committee (Centre Hospitalier Intercommunal de Créteil, France) and was

registered at ClinicalTrials.gov: NCT0441231. Legal guardians and children were informed with a written non-opposition form. Written informed consent was not required to participate in this study in accordance with the national legislation and the institutional requirements.

Author contributions

RC, SB, CJ, and CL designed the study. CA, AF, FC-S, OR, AA, CB, and BV made the acquisition of the study data. RC, AR, SB, and CL drafted the initial manuscript and agree to be accountable for all aspects of the work in ensuring that questions related to the accuracy or integrity of any part of the work are appropriately investigated and resolved. All the authors analyzed the data, revised critically the manuscript for important intellectual content, and provide approval for publication of the content. All authors contributed to the article and approved the submitted version.

Funding

This study has been self-funded by ACTIV, including the purchase of the tests.

Acknowledgments

We are grateful to the investigators of the VIGIL Study Network: Ansoborlo Sophie; Auguste Anna; Auvrignon Anne; Batard Christophe; Beaufils Philippe Florence; Bellemin Kaarine; Benani Mohammed; Bonnel Anne-Sophie; Boulanger Sophie; Burtscher Alain; Cahn Sellem Fabienne; Cambier Nappo Eliane; Chartier Albrech Chantal; Chevé Anne; Cheymol Jacques; Cohen Robert; Come Matthieu; Condor Roxana; Corrard François; Cosson Marie-Anne; Cottias Anne-sylvestre; Coudy Caroline; d'Ovidio Nadia; Dagrenat Véronique; De Pontual Loic; Deberdt Patrice; Defives Isabelle; Delavie Nadège; Delvaux Samantha; Desandes Roxane; Desvignes Veronique; El Aouane El Ghomari Imane; Elbez Annie; Elmerich Florence; Favier Marion; Filleron Anne; Gebhard Francoise; Gelbert Nathalie; Givois Annick; Gorde Stephanie; Grué-Fertin Pascaline; Guiheneuf Cécile; Haas Hervé; Hassid Frédéric; Hennequin Stéphanie; Hubinois Sylvie; Jhaouat Imen; Jouty Cécile; Karaa Daniele; Kherbaoui Louisa; Kochert Fabienne; Koskas Marc; Langlais Sophie; Laporte Eve; Le Mouel Fanny; Lemaître Chloé; Louvel Murielle; Magendie Christine; Mazaud Sabine; Mercier Antoine; Merckx Audrey; Milliard Dominique; Minette Delphine; Mizzi-Rozier Marie; Monteil Stéphane; Nold Bénédicte; Perrier Coline; Pflieger Hugues; Picard Karine; Pressac Isabelle; Ravilly Sophie; Roche Christine; Romain Olivier; Romano Stephane; Roques Gaëlle; Savajols Elodie;

References

1. European Centre for Disease Prevention and Control. *Data on Country Response Measures to COVID-19.* (2022). Available online at: https://www.ecdc. europa.eu/en/publications-data/download-data-response-measures-covid-19 (accessed September 6, 2022).

2. Mercer TR, Salit M. Testing at scale during the COVID-19 pandemic. Nat Rev Genet. (2021) 22:415–26. doi: 10.1038/s41576-021-00360-w

3. Haute Autorité de Santé. *COVID-19 - Quel Test Faire?* (2022). Available online at: https://www.has-sante.fr/upload/docs/application/pdf/2022-01/fu_covid19_ quel_test_faire_19_janvier.pdf (accessed September 6, 2022).

4. Direction de la Recherche, des Etudes, de l'Evaluation et des Statistiques. 168 215 000 tests RT-PCR et antigéniques validés en 2021. (2022). Available from: https://drees.solidarites-sante.gouv.fr/delais-covid19-2022-01-06#:\$\sim\$: text=Cela%20repr%C3%A9sente%201%20419%20700,particuli%C3%A8rement %20les%2016%2D65%20ans (accessed September 6, 2022).

5. Mina MJ, Parker R, Larremore DB. Rethinking Covid-19 test sensitivity - a strategy for containment. *N Engl J Med.* (2020) 383:e120. doi: 10.1056/NEJMp2025631

6. Fujita-Rohwerder N, Beckmann L, Zens Y, Verma A. Diagnostic accuracy of rapid point-of-care tests for diagnosis of current SARS-CoV-2 infections in children: a systematic review and meta-analysis. *BMJ Evid Based Med.* (2022). doi: 10.1101/2021.08.11.21261830

7. Jung C, Levy C, Varon E, Biscardi S, Batard C, Wollner A, et al. Diagnostic accuracy of SARS-CoV-2 antigen detection test in children: a real-life study. *Front Pediatr.* (2021) 9:647274. doi: 10.3389/fped.2021.647274

8. Palmas G, Moriondo M, Trapani S, Ricci S, Calistri E, Pisano L, et al. Nasal swab as preferred clinical specimen for COVID-19 testing in children. *Pediatr Infect Dis J.* (2020) 39:e267–e70. doi: 10.1097/INF.00000000002812

9. Klein JAF, Kruger LJ, Tobian F, Gaeddert M, Lainati F, Schnitzler P, et al. Head-to-head performance comparison of self-collected nasal versus

Sellam Aurelie; Seror Elisa; Somerville David; Stock Claire; Thiebault Georges; Thollot Franck; Tizi Oualou Linda; Touitou Robert; Tourneur Florence; Turberg-Romain Catherine; Vassal Melanie; Vie Le Sage Francois; Vigreux Jean-Christophe; Virey Brigitte; Werner Andreas; Wollner Alain; Zenatti Isabelle; and Zouari Morched, ACTIV team: Borg Marine; Prieur Aurore; Ramay Isabelle, and the CRC team: Brussieux Maxime and Hoffart Cécile from the Clinical Research Center of the CHI Créteil.

Conflict of interest

The authors declare that the research was conducted in the absence of any commercial or financial relationships that could be construed as a potential conflict of interest.

Publisher's note

All claims expressed in this article are solely those of the authors and do not necessarily represent those of their affiliated organizations, or those of the publisher, the editors and the reviewers. Any product that may be evaluated in this article, or claim that may be made by its manufacturer, is not guaranteed or endorsed by the publisher.

professional-collected nasopharyngeal swab for a WHO-listed SARS-CoV-2 antigen-detecting rapid diagnostic test. *Med Microbiol Immunol.* (2021) 210:181–6. doi: 10.1007/s00430-021-00710-9

10. Lindner AK, Nikolai O, Kausch F, Wintel M, Hommes F, Gertler M, et al. Head-to-head comparison of SARS-CoV-2 antigen-detecting rapid test with self-collected nasal swab versus professional-collected nasopharyngeal swab. *Eur Respir J.* (2021) 57:2003961. doi: 10.1183/13993003.03961-2020

11. Lindner AK, Nikolai O, Rohardt C, Kausch F, Wintel M, Gertler M, et al. Diagnostic accuracy and feasibility of patient self-testing with a SARS-CoV-2 antigen-detecting rapid test. *J Clin Virol.* (2021) 141:104874. doi: 10.1016/j.jcv.2021.104874

12. Ricci S, Lodi L, Citera F, Nieddu F, Moriondo M, Guarnieri V, et al. How home anterior self-collected nasal swab simplifies SARS-CoV-2 testing: new surveillance horizons in public health and beyond. *Virol J.* (2021) 18:59. doi: 10.1186/s12985-021-01533-z

13. Center for Disease Control and Prevention. *Interim Guidelines for Collecting and Handling of Clinical Specimens for COVID-19 Testing.* (2022). Available online at: https://www.cdc.gov/coronavirus/2019-ncov/lab/guidelines-clinical-specimens.html (accessed September 6, 2022).

14. Institut Pasteur. Cahier des charges pour le séquençage du SARS-CoV-2. (2022). Available online at: https://www.pasteur.fr/fr/sante-publique/centres-nationauxreference/cnr/virus-infections-respiratoires-dont-grippe.

15. Dalmat Y-M. Autotest Covid-19 AAZ : la preuve par la PCR. Option/Bio. (2021) 32:6-7. doi: 10.1016/S0992-5945(21)00091-X

16. Santé Publique France. Enquêtes Flash : évaluation de la circulation des variants du SARS-CoV-2 en France. (2022). Available online at: https://www.santepubliquefrance.fr/etudes-et-enquetes/enquetes-flash-evaluation-de-lacirculation-des-variants-du-SARS-CoV-2-en-france#block-337272 (accessed March 3, 2022). 17. Lopera TJ, Alzate-Angel JC, Diaz FJ, Rugeles MT, Aguilar-Jimenez W. The usefulness of antigen testing in predicting contagiousness in COVID-19. *Microbiol Spectr.* (2022) 10:e0196221. doi: 10.1128/spectrum.01962-21

18. Jefferson T, Spencer EA, Brassey J, Heneghan C. Viral Cultures for coronavirus disease 2019 infectivity assessment: a systematic review. *Clin Infect Dis.* (2021) 73:e3884–e99. doi: 10.1093/cid/ciaa1764

19. Société Française de Microbiologie. Avis du 25 septembre 2020 de la Société Française de Microbiologie (SFM) relatif à l'interprétation de la valeur de Ct (estimation de la charge virale) obtenue en cas de RT-PCR SARS-CoV-2 positive sur les prélèvements cliniques réalisés à des fins diagnostiques ou de dépistage. (2021). Available online at: https://www.sfm-microbiologie.org/wp-content/ uploads/2021/01/Avis-SFM-valeur-Ct-excre%CC%81tion-virale-_-Version-def-14012021_V4.pdf (accessed September 6, 2022).

20. Tom MR, Mina MJ. To Interpret the SARS-CoV-2 test, consider the cycle threshold value. *Clin Infect Dis.* (2020) 71:2252–4. doi: 10.1093/cid/ ciaa619

21. McCulloch DJ, Kim AE, Wilcox NC, Logue JK, Greninger AL, Englund JA, et al. Comparison of unsupervised home self-collected midnasal swabs with clinician-collected nasopharyngeal swabs for detection of SARS-CoV-2 infection. *JAMA Netw Open.* (2020) 3:e2016382. doi: 10.1001/jamanetworkopen.2020. 16382

22. Gras-Le Guen C, Cohen R, Rozenberg J, Launay E, Levy-Bruhl D, Delacourt C. Reopening schools in the context of increasing COVID-19 community transmission: The French experience. *Arch Pediatr.* (2021) 28:178–85. doi: 10.1016/j.arcped.2021.02.001

23. Center for Disease Control and Prevention. *What to Know About COVID-19 Testing in Schools.* (2022). Available online at: https://www.cdc.gov/coronavirus/2019-ncov/community/schools-childcare/what-you-should-know.html (accessed September 6, 2022).

24. Service-Public.fr. Le protocole sanitaire applicable dans les écoles, collèges et lycées à partir du 14 mars 2022. (2022). Available online at: https://www.service-public.fr/particuliers/actualites/A15074 (accessed September 6, 2022).

25. Soni A, Herbert C, Filippaios A, Broach J, Colubri A, Fahey N, et al. Comparison of rapid antigen tests' performance between Delta (B.1.61.7; AY.X) and Omicron (B.1.1.529; BA1) variants of SARS-CoV-2: secondary analysis from a serial home self-testing study. *medRxiv*. 2022:2022.02.27.22271090. doi: 10.1101/2022.02.27.22271090

26. Gentile A, Juarez MDV, Lucion MF, Pejito MN, Alexay S, Orqueda AS, et al. COVID-19 in Children: correlation between epidemiologic, clinical characteristics, and RT-qPCR cycle threshold values. *Pediatr Infect Dis J.* (2022) 41:666–70. doi: 10.1097/INF.00000000003564